




**LANDRENEAU**  
**THORACOSCOPY INSTRUMENTS BY PILLING®**  
Specialty instruments for video-assisted thoracoscopic procedures

*Standard open hand instrumentation versus the Landreneau Instrument*

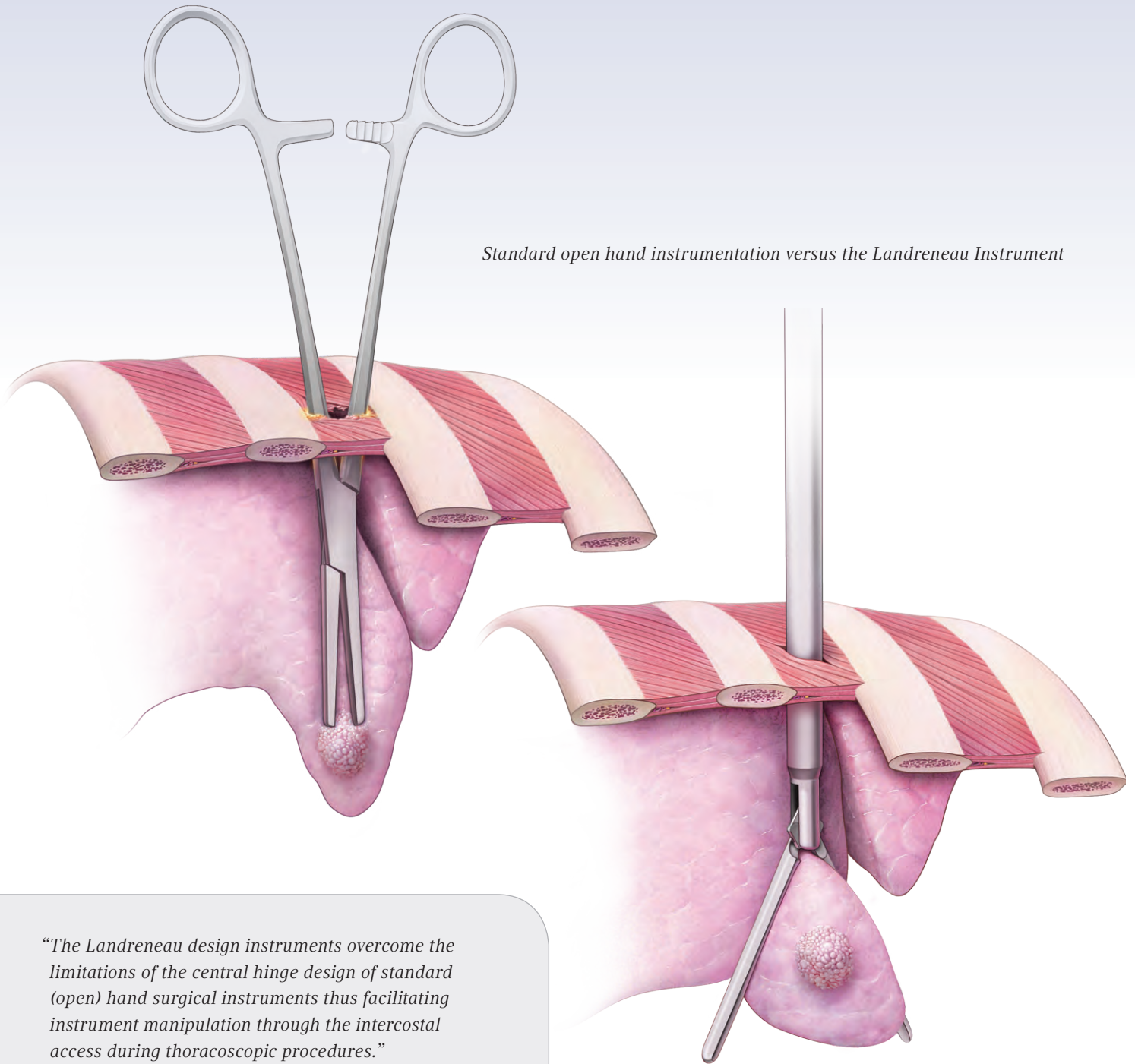

*“The Landreneau design instruments overcome the limitations of the central hinge design of standard (open) hand surgical instruments thus facilitating instrument manipulation through the intercostal access during thoracoscopic procedures.”*

Rodney J. Landreneau, M.D.  
University of Pittsburgh Medical Center

Teleflex Incorporated (NYSE: TFX) is a diversified global company, distinguished by a significant presence in healthcare, with niche businesses that also serve the aerospace and commercial markets.

Teleflex Medical is committed to partnering with health-care providers in critical care and surgery to provide solutions that help reduce infections and improve patient and provider safety. The company also produces surgical instruments and devices, cardiac devices and other specialty products for device manufacturers.

The Teleflex Medical family of brands includes ARROW<sup>®</sup>, BEERE<sup>®</sup>, DEKNATEL<sup>®</sup>, GIBECK<sup>®</sup>, HUDSON RCI<sup>®</sup>, KMEDIC<sup>®</sup>, PILLING<sup>®</sup>, PLEUR-EVAC<sup>®</sup>, RUSCH<sup>®</sup>, SHERIDAN<sup>®</sup>, SMD<sup>™</sup>, TAUT<sup>®</sup>, TFX OEM<sup>®</sup> and WECK<sup>®</sup>, all of which are trademarks or registered trademarks of Teleflex Incorporated or its affiliates.

For detailed information, see [www.teleflexmedical.com](http://www.teleflexmedical.com)

#### **TELEFLEX MEDICAL**

2917 Weck Drive  
Research Triangle Park, NC 27709


#### Mailing Address

Teleflex Medical  
PO Box 12600  
Research Triangle Park, NC 27709  
Toll Free 866.246.6990  
Phone 919.544.8000  
Intl 919.361.3961



211500  
LANDRENEAU THORACOSCOPY TRAY

211600  
LANDRENEAU THORACOSCOPY INSTRUMENT SET

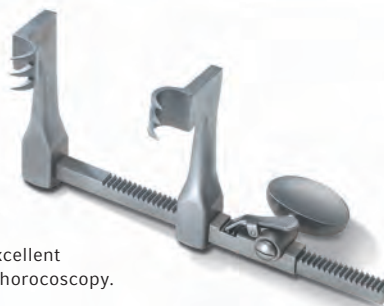

341310  
HEITMILLER MINI  
THORACOTOMY  
RETRACTOR  
Optional instrument  
for the Landreneau  
Thoracoscopy  
Instrument Set

Small "footprint" provides excellent exposure for video-assisted thoracoscopy.

## Landreneau Thoracoscopy Instruments

are used in video-assisted thoracoscopy and minimally invasive procedures. Their unique design overcomes the limitations of the central-hinge on standard open instruments, enabling the instruments to be fully functional when placed through a small incision. Elongated shafts and large jaws provide extended reach while the ring-handled design offers the comfort and familiarity of a traditional open instrument. A custom sterilization tray is available to protect your instrument investment.



131114A  
PENNINGTON  
GRASPING FORCEPS  
10mm Diameter  
13cm shaft length,  
locking



211116A  
STRAIGHT MASHER  
GRASPING FORCEPS  
Thin Teardrop,  
10mm Diameter  
21cm shaft length,  
locking



211117A  
CURVED MASHER  
GRASPING FORCEPS  
Thin Teardrop, not  
designed for use with  
a 10mm cannula  
21cm shaft length,  
locking



211118A  
STRAIGHT MASHER  
GRASPING FORCEPS  
Round, not designed  
for use with a 10mm  
cannula  
21cm shaft length,  
locking



211603A  
**STRAIGHT DEBAKEY  
 GRASPING FORCEPS**  
 10mm Diameter  
 21cm shaft length,  
 locking

211245A  
**LAWRENCE FINE  
 DISSECTING FORCEPS**  
 45°curve, not designed  
 for use with a 10mm  
 cannula  
 21cm shaft length,  
 locking

211290A  
**MEEKER DISSECTING  
 FORCEPS**  
 90°curve, not designed  
 for use with a 10mm  
 cannula  
 21cm shaft length,  
 locking

211325A  
**METZENBAUM  
 SCISSORS**  
 Curved, 10mm diameter  
 21cm shaft length,  
 non-locking

215558A  
**STRAIGHT SINGLE  
 ACTION NEEDLE  
 HOLDER**  
 5mm diameter  
 21cm shaft length,  
 locking